FACE MASK ASTM F2100 / ASTM F1862 Fluid Resistance with synthetic blood Instructions for use / 使用說明 EN14683:2019 Annex B Delta P < 6.0mm H₂O Level 2 ASTM F2101 Bacterial Filtration Effciency (BFE) ASTM F2299 Submicron Particle Filtration (PFE) $\geq 99\%$ 16 CFR Part 1610 Class 1 and pre-bend the your face and pull piece to tighten nose piece in the loops behind the fit. the middle. 拉開口罩並預先彎曲 將口罩覆蓋面部及掛 調校鼻樑壓條至合適 Made in Hong Kong 香港製造 鼻樑壓條至中間位置 上耳帶

20 bcs

PRECAUTIONS:

Manufactured by:

mm 26 x 271

з Гауегѕ

S level 00127 MTSA

注意事項:

(1) This product is a disposable mask (1) 本產品是一次性使用口罩, Please do not reuse, wash or modify 請不要重覆使用、洗滌或改裝 , (2) 建議佩戴連續4小時後 (2) Not recommended to use

continuously for more than 4 hours 需要更換 (3) 請不要在火邊使用 (4) 請勿將本產品存放在潮濕的地方

(3) Do not use it near fire (4) Keep in a cool and dry place (5) 請將本產品存放在兒童接觸不到的地方

(5) Keep out of the reach of children (6) Dispose mask properly after used

Sun Shine Investment Holdings Limited

Unit B, 2/F., Chaiwan Industrial Centre

Email: info@vheart.com.hk

20 Lee Chung Street Chaiwan Hong Kong

(6) 使用後的口罩,請妥善處理 (7) Not made with latex (8) Shelf life 3 years

(7) 本品不含橡膠 (8) 保質期3年

Expiry date / 到期日:

ASTM F2100 Level 2


3 Layers

175 x 95 mm





50 pcs

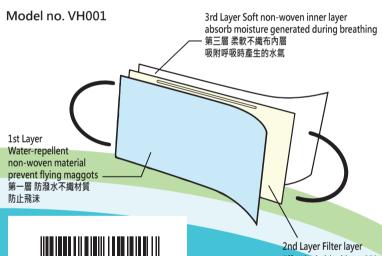

2nd Layer Filter layer Effectively blocking ≥99% of bacteria and particles 第二層 過濾層